



CRYOSPACE" **INNOVATIVE SOLUTIONS** AIR CIRCULATION SYSTEM PATENT PENDING

MEDICAL CERTIFICATE



CRYOSPACE MED

APPLIES QUALITY MANAGEMENT IN DESIGN, PRODUCTION, DISTRIBUTION AND SERVICING OF WHOLE-BODY CRYOTHERAPY CHAMBERS


> COMPLIES WITH THE REQUIREMENTS CONCERNING MEDICAL DEVICES SPECIFIED IN FOLLOWING NORMS: EN ISO 13485:2012 EN ISO 13485:2012/AC:2012





UNIVERSAL MODEL WITH AUTOMATIC ELECTRIC LIFT



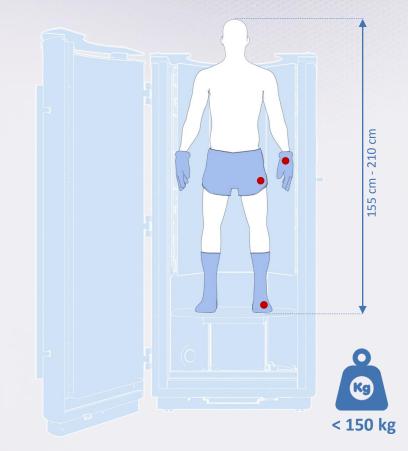

ADJUSTED TO A WIDE VARIETY OF CUSTOMERS


CUSTOMERS HEIGHT from **155 cm** to **210 cm**

CUSTOMERS WEIGHT up to **150 kg**

> SAFETY CLOTHES obligatory use of:


•cotton underwear •gloves •socks •slippers



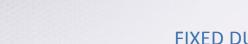




QUALIFIED STAFF NEEDS TO BE PRESENT DURING SESSIONS

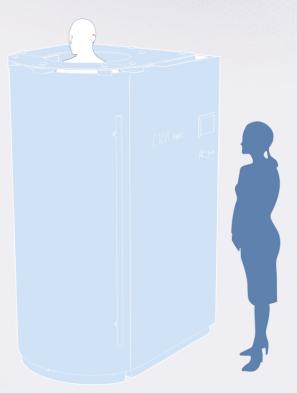

FIXED DURATION OF THERAPY SESSIONS



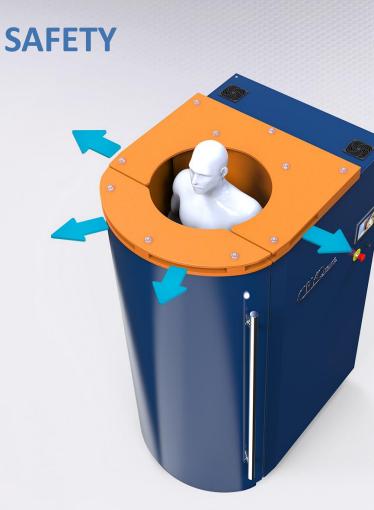



SPECIAL COLLAR DESIGN DRAINS NITROGEN VAPOUR AWAY FROM THE CLIENT





TEMPERATURE STABILITY THANKS TO THE SYSTEM OF AIR CIRCULATION **NITROGEN SUPPLY** TO THE CABIN IN OPTIMAL STATE AND TEMPERATURE **AUTOMATIC EXTRACTION** OF EVAPORATED NITROGEN AFTER SESSION





A DISSASSEMBLED APPLIANCE WILL FIT THROUGH A 80 CM WIDE DOORWAY

PERFORMANCE



	CRYOSPACE	OTHER CRYOGENIC DEVICES
AVERAGE NITROGEN CONSUMPTION	~1kg/1MIN	1,5kg – 3 kg/1MIN
AVERAGE ENERGY CONSUMPTION	0,382 KWh	1,2 – 1,5 kWh
COEFFICIENT OF HEAT PENETRATION FOR	U=0,024 W/(m ^{2.} K)	U=0,035-0,045 W/(m ^{2.} K)
DEVICE SHELL		



PERFORMANCE



CRYO SHEE

CERTIFIED MATERIALS DESIGNATED FOR CRYOGENIC TEMPERATURES

APPROPRIATE AND CONSTANT TEMPERATURE DURING SESSION

-140°C





CRYO_{SPACE}.

CRYOseace



cryospace@cryospace.eu

www.cryospace.eu